

45 Liter Insufflator

Finding Pneumo

Pneumo Sure™



Outstanding Performance

Market Leading Technology

Pneumo Sure

Stryker's 45L Pneumo Sure provides outstanding performance, real time pressure sensing and multiple user modes. With its new integrated tube design, the Pneumo Sure is designed for management of optimal pressure and flow rates to achieve greater accuracy during procedures. The 45L insufflator is designed to handle a wide variety of cases including pediatrics, bariatrics, and vein harvesting. Its ability to safely maintain pneumoperitenium in a range of conditions brings insufflation in the OR to a new level.




Pneumo Sure provides a wide variety of features:

- Multiple User Modes
- Real Time Pressure Sensing
- Software Upgradeable
- Increased Accuracy
- Interactive user Interface

New Touchscreen Interface

Quick & Easy Navigation

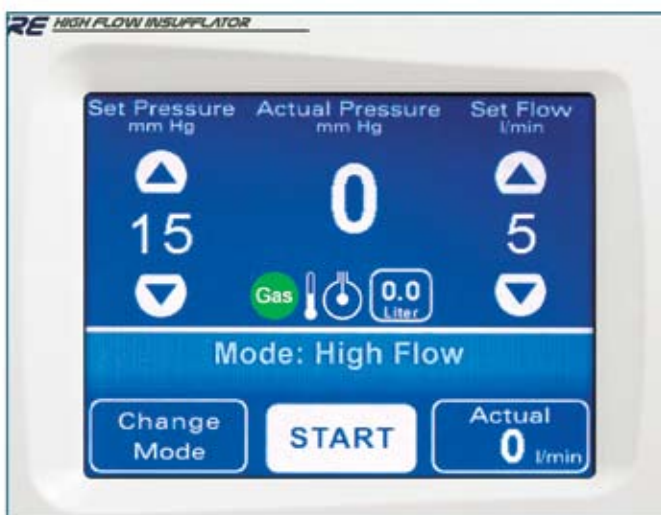

Touchscreen Interface

The Pneumo Sure's intuitive touchscreen ensures quick and easy navigation through the user menu. The redesigned full color interface was specially designed for surgeons and staff.

Features:

- Higher actual flow rate
- Full color touch panel allows for better viewing angles
- SIDNE compatibility allows surgeons to control the insufflator
- Change Modes with the touch of a button

Software Upgrade

The Pneumo Sure was designed to handle the future needs of the OR through an open architecture software platform.

Superior Performance

Across Multiple Specialties

Exclusive Tubeset Design

The Pneumo Sure offers 3 different tubesets with the option of using real time pressure for more accuracy and increased flow. With the advancement of RTP, the insufflator is able to continuously pump CO₂ into the abdomen while constantly measuring pressure.



High Flow Tubing

- Increased Flow
- Integrated Filter
- Ideal for Laparoscopic Cases

High Flow With RTP

- Increased accuracy
- Increased flow
- Ideal for advanced laparoscopic cases



Heated With RTP

- Heated CO₂
- Maintains higher flow rates
- Increased accuracy



System Specifications

45 Liter PNEUMO SURE™	SYSTEM SPECIFICATIONS
Power Supply Input:	100-240 Vac; 3.15A
Frequency:	50-60 Hz
Power Consumption:	150VA
Dimensions:	5.83" (h) x 18.7" (d) x 12.52"(w)
Weight:	19.84 lbs
Operation Conditions - temp:	50-104 F
Operation Conditions - humidity:	30%-75%

Item code	Description
Pneumo Sure High Flow Insufflator (includes Pediatric and High Flow modes)	
620-040-601	Incl. manual: EN,DE,FR,NL,PL Sidne
620-040-602	Incl. manual: EN,NO,SV,DA,FI Sidne
620-040-603	Incl. manual: EN,ES,IT,PT,GR,RO Sidne
Pneumo Sure XL (includes Pediatric, High Flow, Bariatric and Vein Harvesting modes)	
620-040-611	Incl. manual: EN,DE,FR,NL,PL Sidne
620-040-612	Incl. manual: EN,NO,SV,DA,FI Sidne
620-040-613	Incl. manual: EN,ES,IT,PT,GR,RO Sidne
Scope of delivery:	
1000-401-039	Manual: EN,DE,FR,NL,PL version depends on model
1000-401-041	Manual: EN,ES,IT,PT,GR,RO ordered (see above)
1000-401-040	Manual: EN,NO,SV,DA,FI
Z0101-01	Power supply cord EU
Z0600-01	Wrench for High-pressure hose US 9/16"
Upgrade Kit:	
620-040-620	Pneumo Sure Mode Upgrade Kit to XL
Tubing sets:	
Z1463-39	Pneumo Sure tube set, disposable (no RTP) (box of 10 pcs.)
Z1461-39	Pneumo Sure tube set, heating with RTP, disposable (box of 10 pcs.)
Z1462-39	Pneumo Sure tube set with RTP, disposable (box of 10 pcs.)
Z1460-39	Pneumo Sure tube set, heating, reusable (no RTP)
Gas connection:	
Z5010-01	Gas connector NIST for central gas supply (EU standard)
High pressure hose for bottle supply:	
Z5044-01	High pressure tube device US bottle DIN
Z5045-01	High pressure tube device US bottle ISO
Z0499-01	High pressure tube device US bottle Pin Index
Z0498-01	High pressure tube device US bottle US
Wrench for High pressure hose:	
Z0600-01	Fork wrench for connection US
Z0601-01	Fork wrench for connection DIN
Z0602-01	Fork wrench for connection Pin Index
Z0603-01	Fork wrench for connection ISO
High pressure hose for central gas supply:	
Z5030-01	High pressure tube for central gas supply 3m, NIST/DIN
Z5031-01	High pressure tube for central gas supply 5m, NIST/DIN
Z5032-01	High pressure tube for central gas supply 3m, NIST/AGA
Z5033-01	High pressure tube for central gas supply 5m, NIST/AGA
Z5034-01	High pressure tube for central gas supply 3m, NIST/British Standard
Z5035-01	High pressure tube for central gas supply 5m, NIST/British Standard
Z5036-01	High pressure tube for central gas supply 3m, NIST/Norme Française
Z5037-01	High pressure tube for central gas supply 5m, NIST/Norme Française
Z5038-01	High pressure tube for central gas supply 3m, NIST/UNI (Italian)
Z5039-01	High pressure tube for central gas supply 5m, NIST/UNI (Italian)
Power Supply Cord:	
Z0101-01	Power cord EURO 2m
Z0176-00	Power cord UK 2m
Z0177-00	Power cord CH 2m
Misc.:	
Z0164-01	Tube connector ISO-F/olive, reusable (spare part for Z1460-39)
Z0701-01	Insufflation needle 100 mm
Z0703-01	Insufflation needle High Flow 100 mm
Z0536-01	Filter ISO-F/ISO-M, Hydrophobic, Disposable (for use with Z1460-39)

stryker[®]

Joint Replacements

Trauma, Extremities & Deformities

Cranio-maxillofacial

Spine

Biologics

Surgical Products

Neuro & ENT

Interventional Pain

Navigation

Endoscopy

Communications

Imaging

Patient Handling Equipment

EMS Equipment

Stryker SA
Cité Centre
Grand-Rue 90
1820 Montreux
Switzerland

t: +41 21 966 12 01
f: +41 21 966 12 00

www.stryker.eu

This document is intended solely for the use of healthcare professionals.

The products listed above are CE marked according to the Medical Device Directive 93/42/EEC.

A surgeon must always rely on his or her own professional clinical judgment when deciding whether to use a particular product when treating a particular patient. Stryker does not dispense medical advice and recommends that surgeons be trained in the use of any particular product before using it in surgery.

The information presented in this brochure is intended to demonstrate a Stryker product. Always refer to the package insert, product label and/or user instructions before using any Stryker product. Products may not be available in all markets.

Product availability is subject to the regulatory or medical practices that govern individual markets. Please contact your Stryker representative if you have questions about the availability of Stryker products in your area.

Stryker Corporation or its divisions or other corporate affiliated entities own, use or have applied for the following trademarks or service marks: Pneumo Sure, SIDNE, Stryker. All other trademarks are trademarks of their respective owners or holders.

Products referenced with [™] designation are trademarks of Stryker.

Products referenced with [®] designation are registered trademarks of Stryker.



Literature Number: MTX1000900581 Rev. A
MTX4075/GS 02/09

Copyright © 2009 Stryker